


The Swann-Morton Surgical Blade Remover

Single use and available sterile in boxes of 50 it can be used within the Operating Room at the point of the procedure to remove the full range of surgical blade shapes from either a No.3 or a No.4 fitment handle. Protecting staff from sharps injuries during contaminated blade removal it is also transparent allowing for the count at the end of the procedure before disposing in an appropriate sharps container.

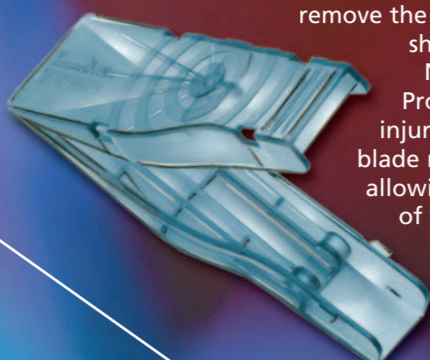


Skin Graft Blade (SG3)

In stainless steel, with 5cm long cutting edge. Fits a No.3 fitment handle. Sterilised, packed individually, boxed in 20's.



The "KLEEN" Blade Management System

Available in a range of popular blade sizes and is compatible with our No.3 and No.4 surgical handles. The cartridge allows for safe attachment and removal of the blade whilst the sheath can be retracted fully exposing the blades during use.

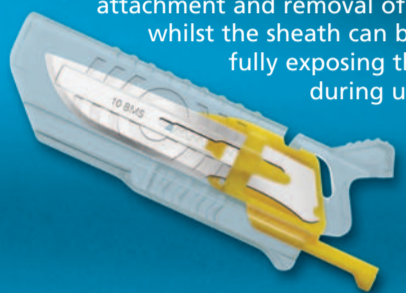

For more information about Cygnetic, KLEEN and our range of Sharps Safety products refer to our website at: www.swann-morton.com
Our suite of training films are available at: <http://www.swann-morton.com/interactive.php>

Cygnetic® (Showing the Attach/Detach position)

Cervical Biopsy Blade
In stainless steel, sterilised, supplied in single plastic trays in boxes of 10. Fits the Long and Short Swann Major Handles which can be viewed on the website along with the Swann-Major blades.

"Minor" Disposables
Can be used as a component part of a "Procedure Kit". Small polystyrene handle can be supplied fitted with blade No's 10, 10A, 11, 15 or 15A. Polyethylene blade guard. Non-sterile and bulk-packed.



Scan this code for more product information

Disposable Scalpels
Available with all standard stainless steel blade shapes, complete with CE marked polystyrene handles with a graduated measure (for indication only), individually sealed in a tray designed to assist with safe passing protocols, sterilized and boxed in units of 10.

Guarded Disposable Scalpels
Available with the most popular Stainless Steel blade shapes, and are fitted with a polyethylene protective guard. The sterile scalpels are individually wrapped in a view pack, and boxed in units of 10. Also available in a bulk, non-sterile format. For use in custom procedure trays and kits.

Disposable Scalpels with Retractable Blade
The Retractable range have colour coded sliders to indicate the blade shape, an audible click when engaging the usage, passing and locked positions and ribbed, slip resistant finger grips for a secure hold in a gloved hand. They are available in Stainless Steel blade shapes 10, 10A, 10R, 11P, 12, 14, 15, 15C, 20, 21, 22, 23, 24 and 25A in sterile packs of 10 or 25. There is also a Stitch Cutter in the range. All shapes are also available in non-sterile packs of 500. For use in custom procedure trays and kits.

"Fine" Blades and Handles
The "Fine" range has small, slim, precise blades in Stainless Steel for delicate operations. Six special handles, including the SF13 and SF23 with independently screwed top sections to facilitate blade-changing, take all blades, supplied in individual metal foil packets, available in dispenser units of 25. Handles are supplied individually in plastic wallets.

Standard Surgical Blades and Handles – No.3 Fitment
Surgical Blade No's 6, 9, 10, 10A, 10R, 11, 11P, E/11, Sabre E/11, 12, 12D, 13, 14, 15, 15A, 15C, 15T, Sabre D/15, 16, 40 fit No's 3, 3L, 5B, 7, 9, B3 handles. All handles are available in Stainless Steel, individually wrapped, and boxed in units of 10. The easily gripped B3 is a light, adaptable stainless handle, supplied individually in plastic wallets as are the bulbous 5B Handles.

Standard Surgical Blades and Handles – No.4 Fitment
Surgical Blade No's 18, 19, 20, 21, 22, 22A, 23, 24, 25, 25A, 26, 27, 36 fit No's 4, 4L and 6B handles. All handles are available in Stainless Steel, individually wrapped, and boxed in units of 10. The larger 6B handle is supplied individually in plastic wallets. Standard surgical blades in carbon steel (sterile or non-sterile) or in stainless steel (sterile). Boxed in 100's (sterile blades individually foil-wrapped; non-sterile sealed in foil packs of 5). Non-sterile blades also available individually packed in peel-open packets, boxed in 100's. Sabre E/11 and D/15 blades are produced in carbon steel only, sterile, boxed in 100's.

Cygnetic® Surgical Blades and Handles
For use within rigorous Orthopaedic and General surgical procedures. Cygnetic® provides the surgeon with additional strength and stability through the patented combination of its specially designed blade and handle. Handles are available in Stainless Steel, Non-Sterile, individually sleeved. Blades shapes 10, 11, 15 and 20 in Stainless Steel, Sterile, individually V.P.I. protected in Metal Foil Packets and boxed in units of 50.

PM40 and PM40B Post Mortem knife and Disposable Scalpels
A Post Mortem knife with a choice of two heavy duty carbon steel blades, the pointed PM40 or bull nosed PM40B. Screw type fitting to a stainless steel handle or complete as a convenient disposable scalpel with a protective guard. Spare blades are individually wrapped and boxed in 10's whilst the scalpels are available in units of 5 and 100.

Post-Mortem PM8 Handle & PM60 and PM60B Blades
Further extending the Post Mortem range, the pointed PM60 and bull nosed PM60B blades offer the ease of attaching to a PM8 handle in a traditional manner. Both blades VPI protected in single-peel packets, boxed in 50's. Handles supplied individually in plastic wallets.

Stitch Cutters
Requiring no handles, the standard stitch cutter is in carbon steel, the midi and long versions in stainless steel. The No 3 stainless stitch cutter can be used in conjunction with a No 3 handle for more aseptic techniques. Available individually foil-wrapped in sterile packs of 100.

Surgical Blades and Handles by Swann-Morton®

www.swann-morton.com



All items are shown actual size



"Swann-Morton" and the "ring pattern" logo are registered trademarks of Swann-Morton Limited and related companies. The products shown are manufactured for use by trained health care professionals and can cause injury if used or handled incorrectly. The information and illustrations used in connection with the products shown on this leaflet are for general guidance only and users must seek proper training and instruction in their use before using or handling such products. Swann-Morton Limited exclude all liability for injury to the extent permitted by law.